## UK – System 75 Wall Mounted ESD

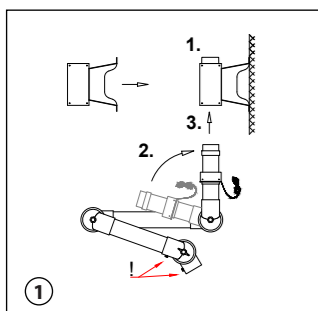
Fig. 1

For mounting in an ESD working area the earthing wire on the socket pipe is very important. Take care not to damage it during the mounting.

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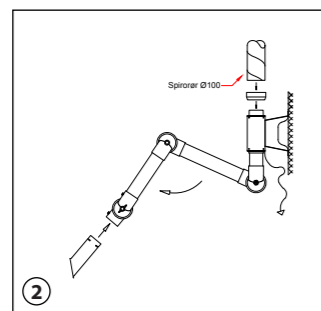
Fig. 2

1. Connect the main ventilation duct directly to the top of the wall bracket. For other connection options, see instruction for wall bracket.

2. In an ESD area a potential equalizing must be ensured for the surrounding working area. The earthing wire pre-attached to the extraction arm must be replaced with the one in the plastic bag marked "ESD". The new earthing wire with a 1MΩ resistance must be connected to the central earthing for the furniture and tools. When turning the extraction arm be aware not to break the 1MΩ resistance on the earthing wire.
3. Unfold the extraction arm in the direction of the arrow.
4. Place the connection pipe of the accessory in the joint. Turn the slit on the connection pipe towards the small fingerscrew on the joint. When screwed in the fingerscrew will hold the accessory in place in the joint.
5. Finally make sure the damper is open.
6. The extraction arm is now ready for use.

## UK – System 75 Ceiling Mounted ESD

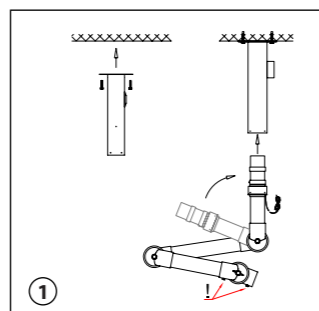
Fig. 1

For mounting in an ESD working area the earthing wire on the socket pipe is very important. Take care not to damage it during the mounting.

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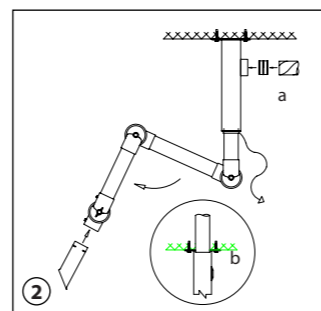
Fig. 2

1. For the side connection (Fig. 2a) the main ventilation duct (80mm dia.) can be connected directly to the outlet on the side of the column. For top connection (Fig. 2b) the main ventilation duct is put through the ceiling and directly into the column.

2. Make sure to tighten the gap between the ceiling and duct with a filler.
2. In an ESD area a potential equalizing must be ensured for the surrounding working area. The earthing wire pre-attached to the extraction arm must be replaced with the one in the plastic bag marked "ESD". The new earthing wire with a 1MΩ resistance must be connected to the central earthing for the furniture and tools. When turning the extraction arm be aware not to break the 1MΩ resistance on the earthing wire. Be aware to break any conductive connection to the main ventilation system.
3. Unfold the extraction arm in the direction of the arrow.
4. Place the connection pipe of the accessory in the joint. Turn the small slit in the connection pipe towards the fingerscrew. When screwed in the fingerscrew will hold the accessory in place in the joint.
5. Finally make sure the damper is open.
6. The extraction arm is now ready for use.

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See MAINTENANCE AND CLEANING at [www.alsident.com](http://www.alsident.com) under "Technical Information"

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## UK – System 75 Wall Mounted EX

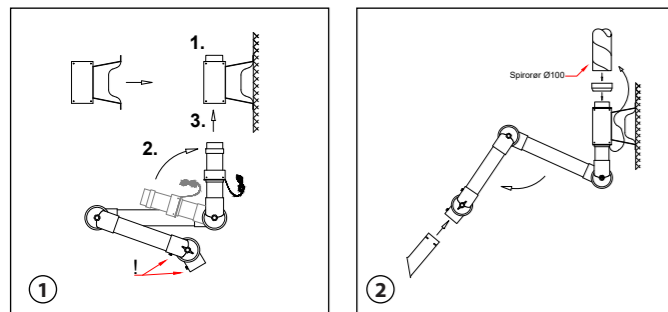
Fig. 1

For mounting in an EX working area the earthing wire on the socket pipe is very important. Take care not to damage it during the mounting.

1. Fasten the wall bracket solidly to the wall or column. When mounting the wall bracket be aware of the considerable strain the operation of the arm will transmit to the bracket.
2. If the wall or column is made of a conductive material the connection to the wall bracket must be non-conductive.
3. Considering the internal spring, it is important to turn the socket pipe in the direction of the arrow (2) and that the fingerscrews (!) face forward.
4. Place the socket pipe in the wall bracket (3) and fasten the square-guide with the included screws.
5. The extraction arm can be turned 190° in the wall bracket.

Fig. 2

1. Connect the main ventilation duct directly to the top of the wall bracket. For other connection options, see instruction for wall bracket.
2. In an explosive working area a potential equalizing must be ensured for the ventilation system. It is therefore important to bond all the conductive parts of the system and to earth them. The extraction arm must be connected to the ventilation duct with the earthing wire pre-attached to the extraction arm. Alsident® System A/S recommends the earthing of the extraction arm to be carried out by personnel experienced in EX areas.
3. Unfold the extraction arm in the direction of the arrow.
4. Place the connection pipe of the accessory in the joint. Turn the slit on the connection pipe towards the small fingerscrew on the joint. When screwed in the fingerscrew will hold the accessory in place in the joint.
5. Finally make sure the damper is open.
6. The extraction arm is now ready for use.

## UK – System 75 Ceiling Mounted EX

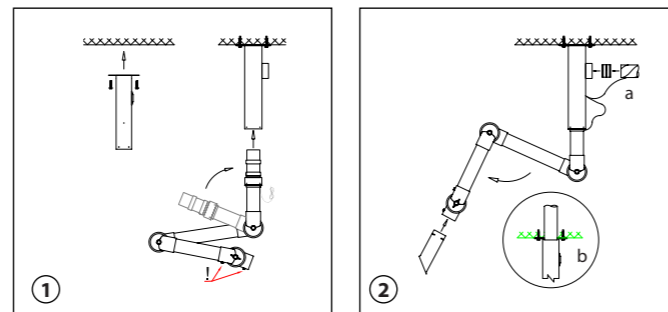
Fig. 1

For mounting in an EX working area the earthing wire on the socket pipe is very important. Take care not to damage it during the mounting.

1. Mount the column directly under the ceiling.
2. The operation of the extraction arm will transmit a considerable strain to the column. The column must therefore be mounted on a solid ceiling.
3. Considering the internal spring, it is important to turn the socket pipe in the direction of the arrow and that the fingerscrews (!) are faced forward.
4. Place the socket pipe in the column and fasten the square-guide with the included screws.
5. The extraction arm can now be turned 360° in the column.

Fig. 2

1. For the side connection (Fig. 2a) the main ventilation duct (80mm dia.) can be connected directly to the outlet on the side of the column. For top connection (Fig. 2b) the main ventilation duct is put through the ceiling and directly into the column. Make sure to tighten the gap between the ceiling and duct with a filler.
2. In an explosive working area a potential equalizing must be ensured for the ventilation system. It is therefore important to bond all the conductive parts of the system and to earth them. The extraction arm must be connected to the ceiling column with the earthing wire pre-attached to the extraction arm. The installer of the ventilation duct must make sure to connect the ceiling column to the ventilations duct with an earthing wire. Alsident® System A/S recommends the earthing of the extraction arm to be carried out by personnel experienced in EX-areas.
3. Unfold the extraction arm in the direction of the arrow.
4. Place the connection pipe of the accessory in the joint. Turn the small slit in the connection pipe towards the fingerscrew. When screwed in the fingerscrew will hold the accessory in place in the joint.
5. Finally make sure the damper is open.
6. The extraction arm is now ready for use.

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## EX

UK: Approved for use in explosive areas and marked

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## ESD

UK:

Approved according to standard IEC 61340-5-1:2016.

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All rights reserved